

BILERUSH SELECTIVE

Biliary Stent for Endoscopy “BILERUSH SELECTIVE”



1. The Smallest Outer Diameter of Delivery System in the world

Developed The Smallest Outer Diameter of Delivery System 5.7Fr(1.9mm)

To use the high rigidity material even though delivery system is the finest delivery system, BSL improved pushability and kink resistance. Physicians can feel the feeling of releasing stent directly.

2. Open Cell Laser Cut Stent


The cell design as “Open Cell” is the widest in the stent which are sold in Japan. The “Open cell” is designed for using multi stenting and Re-intervention

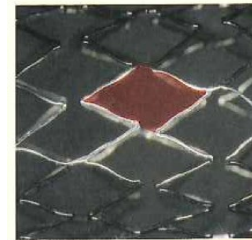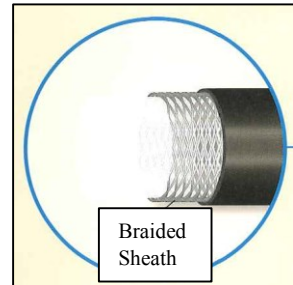


Stent Design



Stent Design of Center part

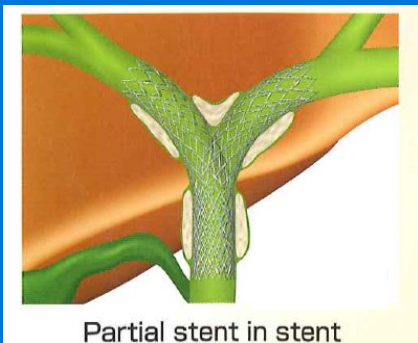

Partial stent in stent

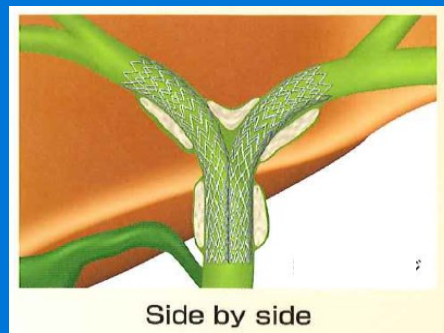


Side by side

3. The expansive force of Stent at beginning of Placement even though AF is low.

The design of open cell to exert the specification of Stent, we have selected the Ni-Ti Alloy which has the stabled expansive force after placement.

Product Feature

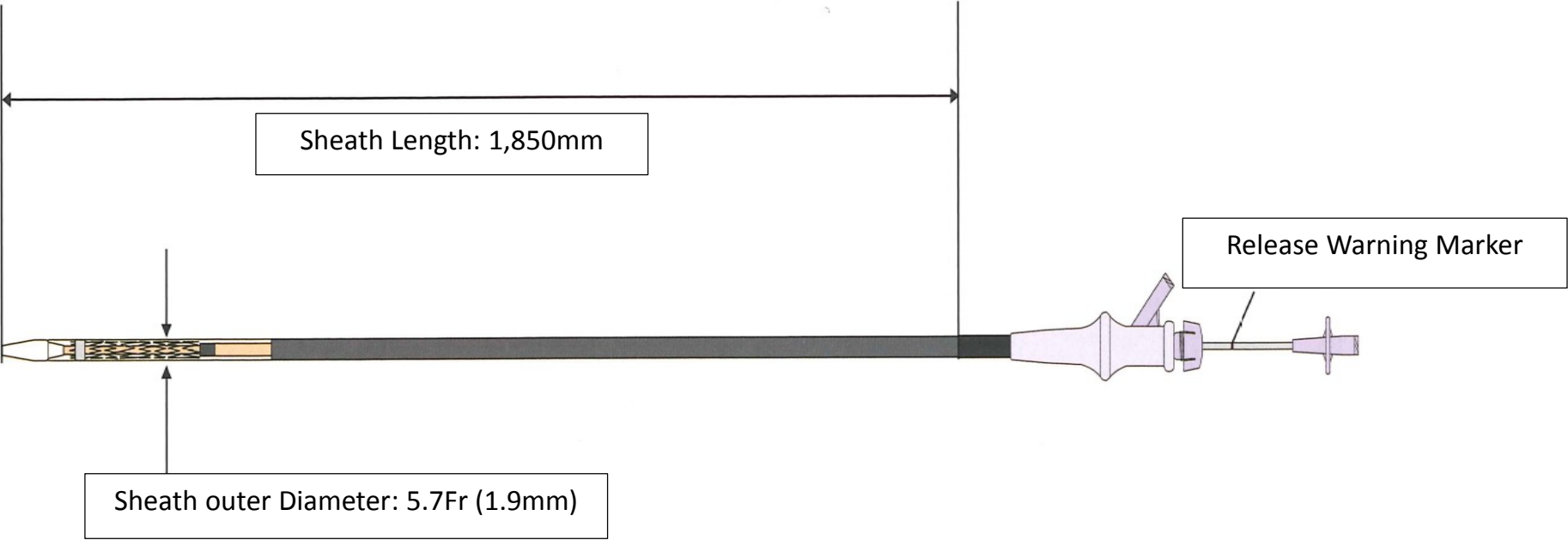

Catalog No.	Stent Diameter	Stent Length	Sheath outer diameter	Guide wire used with	Sheath length
BSL08-40	8mm	40mm	5.7Fr(1.9mm)	0.035inch/0.89mm	1,850mm
BSL08-60		60mm			
BSL08-80		80mm			
BSL10-40	10mm	40mm			
BSL10-60		60mm			
BSL10-80		80mm			