



Product Specification

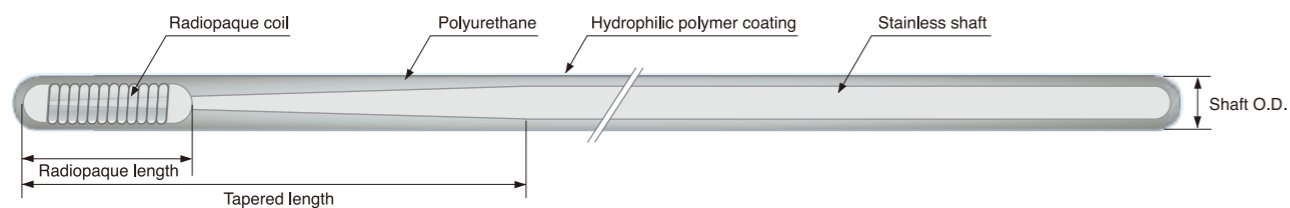

Catalogue Number	Total Length	Tip Shape	Shaft O.D.	Radiopaque Length	Tapered Length
RUN135-16	135cm	Straight	0,41mm (0.016inch)	3cm	40cm
RUN165-16	165cm				
RUN165-16A		Angle(d)			
RUN180-16	180cm	Straight			
RUN180-16A		Angle(d)			

- Ethylene oxide gas sterilization
- Contents : Guidewire, Shaper, Torque Device, Inserter
- For single use only and do not reuse
- Consult instructions for use

Tip Image

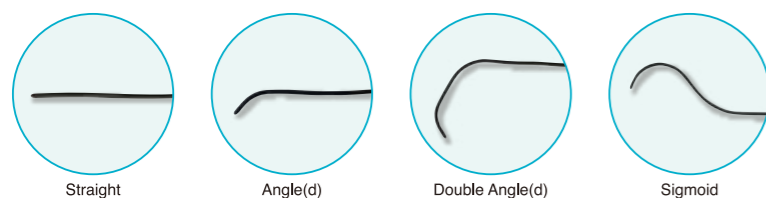

High Shapability / High Re-shapability

Made in Japan

Run&Run

Micro Guidewire

Dedicated for Targeted Chemoembolization





PIOLAX MEDICAL DEVICES, INC.

[Overseas Operations Group]
2F Pacific Marks Shin-Yokohama, 3-18-20 Shin-Yokohama,
Kouhoku-ku, Yokohama-shi, Kanagawa 222-0033 Japan
Tel : +81 45 548 9888 / Fax : +81 45 476 5757
<http://www.piolax-md.co.jp/en/>

"Reshape" is a key for successful procedures.

Shapability

- Same technology in Japanese swords.
- High Shapability by "Ribbon processing".
- High Re-shapability.

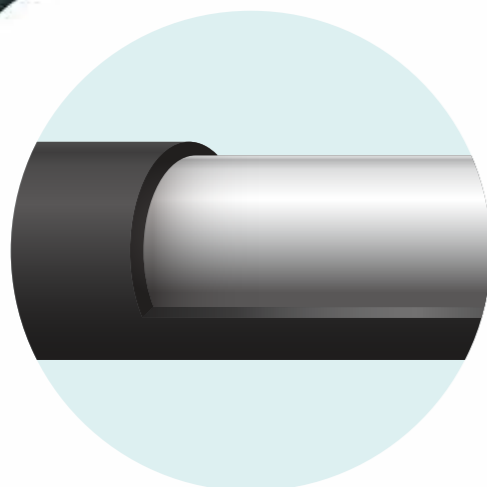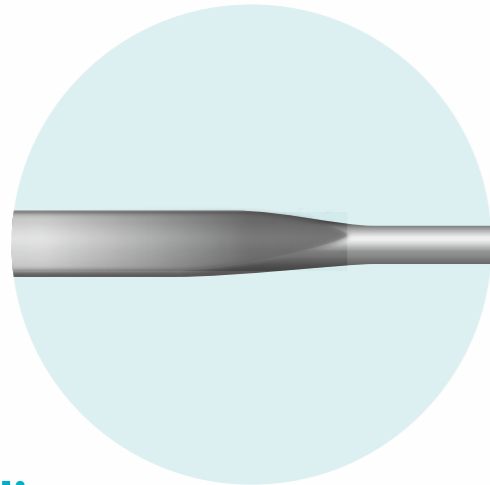


Stainless Core

- Appropriate stiffness, excellent supportability and torquability.

Optimal Clearance

- Leads to ideal trackability of micro catheters.

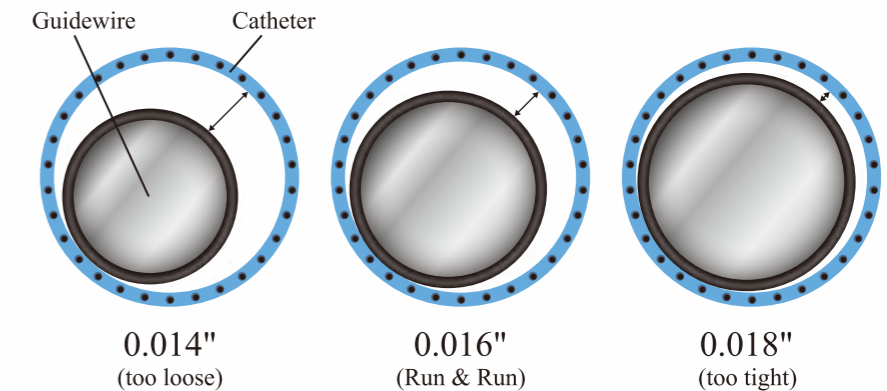


Lubricity

Full coating of Hydrophilic polymer.

Visibility

Reduces intraoperative stress.

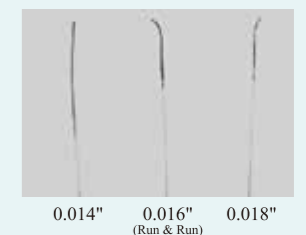

Stability

Decrease the need of additional guidewires.